DE Ferninfrarot Heizkissen  
GB Far Infrared Heating Pad

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**Instruction Manual**

**Modell: HK610**

CE

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Deutschland

## Contents

1.	Important instructions.....	6
2.	Included in delivery .....	7
3.	Product features.....	7
4.	Operating.....	7
5.	Cleaning and maintenance.....	8
6.	Disposal .....	8
7.	Technical data.....	9
8.	Warranty.....	9

### 1. Important instructions

- **Failure to observe the following information can cause personal injury or property damage.**
- **Please DO NOT use the heating pad while sleeping!**
- **Do not use the heating pad on children, disabled or sleeping people, or people who are insensitive to heat and who cannot react to overheating.**
- Children are not allowed to play with the heat pad.
- This heating pad is not intended for use in hospitals.
- **Do not use folded or pushed together.**
- Do not pull, twist or sharply kink the cable.
- **Only use the heat pack with the power adapter supplied.**
- If the heat pad is used for several hours, we recommend setting the lowest temperature level to avoid overheating of the heated body part.
- **The heating pad switches off automatically after 120 minutes.**
- **Discontinue use immediately if the belt becomes uncomfortably hot.**


## 2 . Package Contents

- 1 x heating pad
- 1 x charger
- 1 x Power Adapter
- 1 x instruction manual
- 1 x wash bag

## 3. Product features

The **Thermrup Far Infrared Heating Pad** provides deep heat for the body. It has 6 adjustable temperature levels. The temperature of the heating pad can be adjusted individually. The special rubber heating element has a large heating surface and provides an even heat distribution. The heating pad can be cleaned in a washing machine.

## 4. Operation

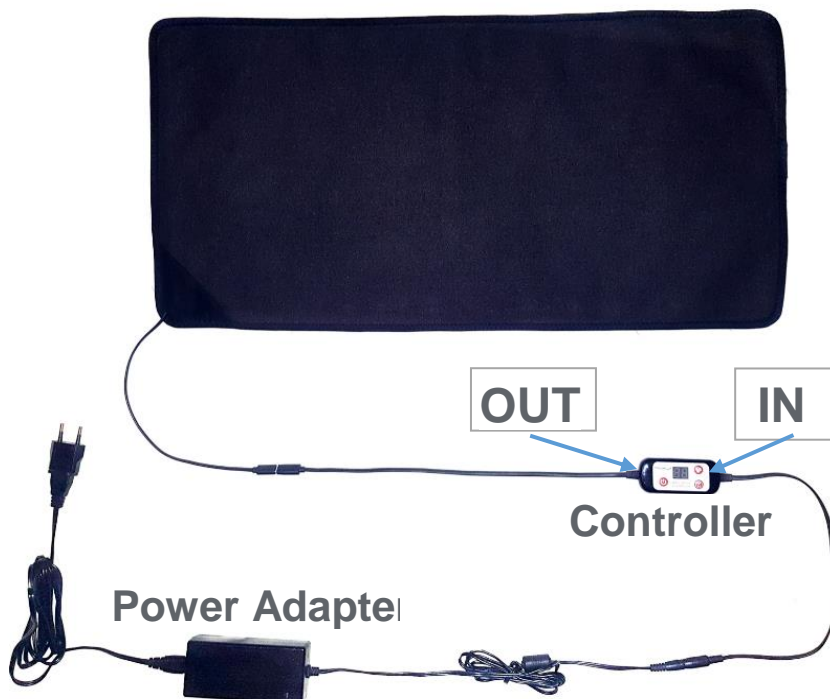- Connect the supplied power supply to the "IN" port of the temperature controller. Next, connect the heating pad to the "OUT" port of the controller.

- Press the power button on the controller to switch on the heating function.
- Then, press the “Temp” button briefly (several times if necessary) to set the required temperature (40-45-50-60-65 ° C).

- Press the power button for approx. 2 seconds to switch off the heating function.
- After use, remove the heating pad from the power supply.
- The heating pad switches itself off automatically after 120 minutes.

**\*\*\* If the heating pad is used for several hours, we recommend setting the lowest temperature level to avoid overheating of the heated body part.**

**\*\*\* Stop using if you see redness on your skin. Excessive heat can cause burns.**

## 5. Cleaning and maintenance

- Before cleaning, disconnect the heating pad from the controller and power adapter.
- Never use aggressive cleaning agents or strong brushes.
- The heating pad can be washed in the washing machine up to 45 ° C in a washing bag or supplied washing bag.
- The heating pad should not be dried or ironed in a dryer.

## 6. Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

## 7. Technical data

Model : HK610  
Power supply  
Heating pad : 12V ; Max. 42W  
Power supply Input : 100-240V ~ 50 / 60Hz  
Output : 12V; 3500mA  
Dimensions : approx. 60x30 cm  
Weight : approx. 415g

## 8. Warranty

We provide a 2-year warranty from the date of purchase for material and manufacturing defects of the product.

The warranty does not apply:

- In the case of damage caused by improper use.
- For wear parts.
- Customer's own fault.
- Once the unit has been opened by an unauthorized workshop.

For the assertion of a warranty within the guarantee period, the customer must prove the purchase.

In case of complaints please contact our service under the following contact:

E-Mail: [service@zewox.de](mailto:service@zewox.de)

### **Zewox GmbH**

Oststr.159

47057 Duisburg

Germany