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www.thermrup.de

Contents

1.	Important instructions.....	6
2.	Included in delivery	7
3.	Product features.....	8
4.	Operating.....	8
5.	Cleaning and maintenance.....	10
6.	Disposal	10
7.	Warranty.....	11

1. Important instructions

- **Non-observance of the following instructions may result in personal injury or damage to property.**
- The following heat belt is intended only for heating the human upper body.
- **Do not use the device on children, disabled or sleeping persons, or heat-insensitive persons who cannot react to overheating.**
- Children shall not play with the appliance.
- **Please do not pull, twist or bend the power lines sharply.**
- The following heat belt is not intended for use in hospitals.
- **When using the heat belt for several hours, we recommend setting the lowest temperature level to avoid overheating the heated part of the body.**
- **Stop using the belt if it becomes uncomfortably hot.**

2. Included in delivery

Please first check the completeness of the device and inspect the product for damage. If the device is damaged, do not operate it and contact your dealer or the service department.

The scope of delivery includes:

- 1 x heated Belt
- 1 x USB charger
- 1 x battery bag
- 1 x Extension USB cable
- 1 x Extension for belt
- 1 x Instruction Manual

3. Product features

The **Thermrup Far infrared Wireless Heated Belt** is designed for indoor and outdoor use provides deep heat with far infrared for back and abdomen with good freedom of movement. The heat belt is easy to put on, very comfortable and safe. It has 3 adjustable temperature settings. The special rubber heating element has large heating area and even heat distribution. The belt can be washed in a machine.


4. Operation

4.1 Operating with power bank

- Now place the heat belt around your hips and abdomen. Next close the Velcro fastener of the belt. The heat belt can be worn in any position on the body at the hips and around the abdomen.
- Set a power bank (not included) in the Battery pack holder of the belt.
- Plug the USB connector of the cable into the USB socket of the powerbank. Then turn on the power bank.
- Stick the battery pouch onto the belt if required.
- Briefly press the power button on the controller to turn on the first level.
- Next, if necessary, press the power key again to set the desired temperature level (1-2-3, red-green-blue). Each single key touch increases the temperature by one level.

The 1st level is the warmest level. (red LED indicator lights up).
The 2nd stage is the middle stage. (green LED indicator lights up).
The 3rd stage is the lowest stage. (blue LED indicator lights up).

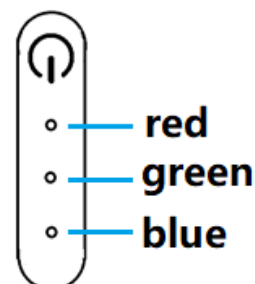- When the blue level is reached and the ON/OFF button is pressed again, the heating function is switched off and the LED indicators no longer light up. After 90 minutes of continuous operation, the heating belt will automatically turn off, restart your powerbank to continue.

!! Discontinue use if you see reddening of the skin or if the belt becomes uncomfortably hot! Excessive heat can cause burns.

4.2 mains operation

- Connect the supplied extension cable to the connection cable of the belt.
- Plug the USB connector of the extension cable into the supplied power supply.
- Press the power button of the controller short to switch on the controller. Then press the power button briefly to set the desired temperature level (1-2-3, red -green-blue).
- Press the power button until all of the LED off to switch off the controller.
- After use, remove the extension cord from the charger.

***** Please use only the supplied charger. *****

5. Cleaning and maintenance

- Before each cleaning, remove the heating pad & the cable from the belt.
- Never use aggressive detergents or strong brushes.
- The heat belt can be washed in the washing machine up to 45°C in a wash bag or in the supplied wash bag.
- The heat belt must not be dried in a dryer or ironed. It should not be used again until the device is completely dry.

6. Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

7. Warranty

We provide a 2-year warranty from the date of purchase for material and manufacturing defects of the product.

The warranty does not apply:

- In the case of damage caused by improper use.
- For wear parts.
- Customer's own fault.
- Once the unit has been opened by an unauthorized workshop.

For the assertion of a warranty within the guarantee period, the customer must prove the purchase.

In case of complaints please contact our service under the following contact:

E-Mail: service@zewox.de

Zewox GmbH

Oststr.159

47057 Duisburg

Germany

www.thermrup.com