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Deutschland

## Contents

1.	Important instructions.....	7
2.	Included in delivery .....	8
3.	Product features.....	8
4.	Operating.....	8
5.	Cleaning and maintenance.....	10
6.	Disposal .....	10
7.	Replacement parts.....	10
8.	Warranty.....	11

### 1. Important instructions

- Non-observance of the following notes may result in personal injury or material damage.
- This heated belt is only designed to warm up the human body.
- **Do not use the heat belt on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).**
- Children shall not play with the appliance.
- The heated belt may only be operated with the battery control unit and the temperature/time adjuster supplied.
- Do not pull, twist or make sharp bends in the cables.
- To charge the battery pack, you must use only the provided charger unit.
- If the heating belt is being used over several hours, we recommend that you set the lowest temperature setting on the battery to avoid overheating the heated body part, which may lead to burns to the skin.
- **Stop using the belt if it becomes uncomfortably hot.**

## 2. Included in delivery

- 1 x heated belt
- 1 x battery
- 1x battery bag
- 1 x charger
- 1 x Extension for belt
- 1 x Instruction Manual

## 3. Product features

The Thermrup **Far infrared Wireless Heated Belt** is designed for indoor and outdoor use provides deep heat with far infrared for back and abdomen with good mobility. It is easy to wear, comfortable and safe. It has 4 adjustable temperature settings. The special rubber heating element has large heating area and even heat distribution. The belt can be washed in washing machine.

## 4. Operation


### 4.1 Charging the battery

- The battery needs to be charged for at least 3 hours prior to initial use or after not having been used for a long time.
- To recharge the battery, connect the cable from the charger to the battery and a power outlet. During charging, the LED will flash. The LED indicator appears constant with the corresponding charging capacity. When the

battery is fully charged, all four LED indicators appear constant. You can then disconnect the charger from the power supply and the battery.

## 4.2 Operating the heated belt

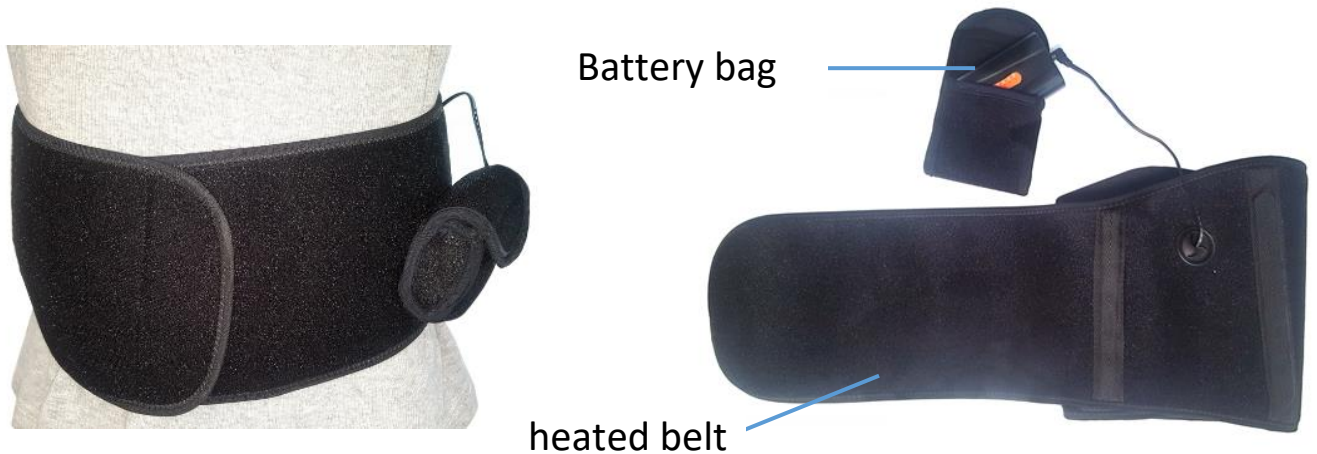- Place the heating belt around your back and abdomen and secure the belt by fastening the Velcro. Depending on the position you put the belt on, it warms different areas of your body
- Connect the power plug of the heated belt and the battery.
- Press the powerbutton for about 2 seconds to turn on the battery. Then press the power button short (if necessary, multiple) to set the desired temperature level (1 to 4). Each press increases the temperature level by one step.
- Puts the battery in the battery bag. Place the battery bag on the heated belt.
- Press the on / off button for about 2 seconds to switch off the battery.

\*\*\*For the application directly with the power supply you need a Temp./Time controller (not included).

**!! Discontinue use if you see reddening of the skin or if the belt becomes uncomfortably hot! Excessive heat can cause burns.**

### 4.3 Checking the battery capacity

Take the cable from the battery and press the power button for 2 second, then the LED will show the remaining capacity of the battery. The LED will light for 3 seconds.

## 5. Cleaning and maintenance

- Before cleaning, take the heating pad from the belt!
- Never use aggressive detergents or strong brushes.
- The heated belt can be washed in the washing machine at 45°C.
- The heated belt must not be dried in a dryer or ironed.

## 6. Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

## 7. Replacement parts

You will find the following accessories and spare parts in our online shop:  
[www.thermrup.com](http://www.thermrup.com)

			
Battery 2200mA	Ersatzakku 4400mA	Extension 20cm	Temperature/time adjuster

## **8. Warranty**

We provide a 2-year warranty from the date of purchase for material and manufacturing defects of the product.

The warranty does not apply:

- In the case of damage caused by improper use.
- For wear parts.
- Customer's own fault.
- Once the unit has been opened by an unauthorized workshop.

For the assertion of a warranty within the guarantee period, the customer must prove the purchase.

In case of complaints please contact our service under the following contact:

E-Mail: [service@zewox.de](mailto:service@zewox.de)

### **Zewox GmbH**

Oststr.159

47057 Duisburg

Germany